




MANAGEMENT SOFTWARE

Optility™

Nuclear Medicine Management Software

Seamless Integration with Capintec™ Dose Calibrators

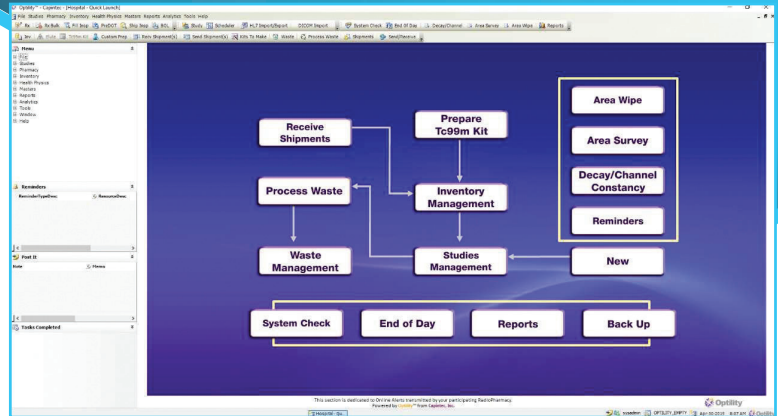

FEATURES

Optility Nuclear Medicine Management Software

- Interface with the full line of Capintec dose calibrators
- HL-7 (HIS) Interface
- DICOM (RIS) Interface
- Bar Coded Inventory
- Online Ordering Functions
- Intuitive User-Friendly Screen Flow
- Health Physics Services
- Highly Customizable Labels and Reports
- HIPAA Compliant with Multiple Security Levels
- Flexible Customer Support Options

Optility is Available in Three Forms:

- **Optility™ NMM** - For hospitals and clinics that dose patients and perform studies.
- **Optility™ CPM** - For radiopharmacies that dispense nuclear or PET prescriptions.
- **Optility™ NPM** - For institutions, such as large universities, that perform studies and dispense prescriptions.

DESCRIPTION


The Optility system is a state-of-the-art Nuclear Medicine Management Software package available in three versions: NMM for Nuclear Medicine imaging facilities that dose patients and perform studies, CPM for commercial pharmacies that dispense prescriptions, and NPM for facilities that perform both operations. It keeps track of all patients, procedures, drugs, inventory and health physics functions with the touch of a few keys.

It allows you to print predefined reports or export the data to Excel spreadsheets or XML formats, and to create your own custom reports.

Protect the important data you've been acquiring for years in the Nuclear Medicine Manager®, Commercial Pharmacy Manager®, or Nuclear Pharmacy Manager® software platform!

With Optility software you can convert all of your information: drugs, procedure settings, health physics configurations, referring physicians, current and historical patient data for immediate use within the program.

OPTIONAL COMPONENTS

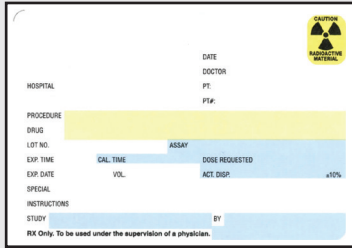

Optility Labels



Dymo Label
Writer Printer 450



Symbol Bar Code
Reader LS 2208



CRC®-55tR Dose
Calibrator

SPECIFICATIONS

Optility Software

Functions

- HL-7 Import and Export Functions
- DICOM Import and Export Functions
- Online Ordering via Web Support Alternate Management Systems
- Export data to Excel or .xlm Formats
- Multiple Security Levels Assure Full Confidentiality of all Data
- Online Software Training

Nuclear Medicine Functions

- Tracks Doses from Orders to Disposition
- Maintains Patient Dose and Exam Records
- Flexible Patient Scheduling Options
- Full Health Physics Compliance
- User-definable Trigger and Notification Levels
- Statistics for Patient Data and Procedures
- Daily Task Log with Action Reminder

Pharmacy Functions

- Tracks and Manages Receiving, Inventory, Components, Kits, and Waste Functions
- Full Traceability from Generator to Unit Dose
- Maintains Batch or Lot Records
- Checks License Limits against Orders
- Business Module includes Invoicing and Financial Reporting
- Analytics to Monitor Productivity and Identify Customer Trends

Services

- Support Knowledge
- Data Imports into Optility Software
- Data Exports out from Optility Software
- Optility Training
- HIS (HL-7) Interface Configuration
- RIS (DICOM) Interface Configuration
- Data migration
- Customization (Labels/Reports)
- Requirements and Analysis
- Software Configuration/enhancement/feature

Requirements

- 64-bit, Windows 10, Windows Server 2016
- Client Server Application
- SQL Server 2017 Express Edition

ORDERING INFORMATION

- Optility Labels: Item #7120-00544
- Dymo Label Writer Printer 450: Item #2000-0007
- Symbol Bar Code Reader LS 2208: Item #5010-0028
- CRC-55tR Dose Calibrator: Item #5130-3234

